

TORQ CONTROL®



FR - CLÉ DYNAMOMÉTRIQUE

EN - TORQUE WRENCH

US - TORQUE WRENCH

DE - DREHMOMENTSCHLÜSSEL

ES - LLAVE DINAMOMÉTRICA

IT - CHIAVE DINAMOMETRICA

PT - CHAVE DINAMOMÉTRICA

NL - TORKSLEUTEL

ZH - 扭矩扳手

AR - مفتاح قياس عزم الدوران

FR - NOTICE D'INSTRUCTION

EN - INSTRUCTIONS FOR USE

US - INSTRUCTIONS FOR USE

DE - GEBRAUCHSANWEISUNG

ES - INSTRUCCIONES DE USO

IT - ISTRUZIONI PER L'USO

PT - MANUAL TÉCNICO DE UTILIZAÇÃO

NL - GEBRUIKSAANWIJZING

ZH - 使用说明

AR - دليل الاستخدام

Anthogyr
A Straumann Group Brand

TABLE OF CONTENTS

- I. DESCRIPTION OF SYMBOLS USED
- II. INDICATIONS FOR USE
- III. AREA OF APPLICATION
- IV. GENERAL SAFETY INFORMATION
- V. TECHNICAL SPECIFICATIONS
- VI. USING THE TORQUE WRENCH
- VII. HYGIENE AND MAINTENANCE
- VIII. CALIBRATION
- IX. REPAIR
- X. GUARANTEES
- XI. MEDICAL DEVICE REFERENCES
- XII. DISPOSAL OF THE PRODUCT

I. DESCRIPTION OF SYMBOLS USED

DANGER SIGNALS

 **WARNING:**
Important Information!

 Wear gloves

REMARKS

 General indications with no danger to people or objects

 Please read this!

 Visual Examination

II. INDICATIONS FOR USE

The torque wrench is used for dental requirements in the area of dental surgery (enamel and dentine treatment, cutting crowns. etc); any other form of use is forbidden and can cause danger. (1) This is an instrument with a precise system of torque control allowing the screwing of prosthetic parts at a definite torque.

III. AREA OF APPLICATION

For this torque wrench, we applied Community Directive 93/42/EEC, transposed by the French regulations on medical equipment (Official Bulletin - L169 - 12 July 1993).

 In accordance with these provisions, the torque wrench must only be used by someone with experience of dental medicine, and specially knowledge in implantology as well as in prosthetic technique related to dental implant for the application described here above, and in compliance with the provisions in effect in relation to the prevention of accidents at work, protection of employment and the instructions for use.

In accordance with these provisions, the user must:

- only use working instruments without defects,
- observe the principle of correct use,
- protect themselves and the patient or third party against any danger
- prevent all contamination by the product.

IV. GENERAL SAFETY INFORMATION

 Before use, check that the torque wrench is not damaged or that no part of it is missing. Use protective gloves, glasses and a mask.

After insertion of a rotary instrument, check that it is correctly held, by a slight rotary movement.

Do not handle the rotary instrument during operation.

If you press the push button during the operation of the instrument, the mandrel can get out of its position. It is recommended to the practitioner to test the torque wrench before

using it on the patient in order to check the good functioning.

 In the event of visible malfunction or damage, immediately stop using the instrument and inform your approved distributor or the manufacturer.

 In the event of any questions on the equipment, contact either your approved distributor or the manufacturer.

V. TECHNICAL CHARACTERISTICS

Weight (g)	130
Input power	Manual
Tool type according to NF EN ISO 1797-1	Type 1
Diameter of tools (mm) according to NF EN ISO 1797-1	2.35
Head height (mm)	12
Head angulation (°)	100
Maximum total length of the tool (mm) *	34
Precision (N.cm)	+/- 1.5

** Approximate values: In the event of the use of rotational instruments which are longer in diameter, the user is responsible for the appropriate choice of operating conditions which will prevent any risk to the patient or to another person.*

Table of the torque value according to Anthogyr prosthetic parts *

In general :

	10	15	20	25	30	32	35
M2 covering screw	✓						
M2 healing screw	✓	✓	✓				
M2 prosthetic screw (gold, titanium)					✓	✓	✓
M2 gingival abutment (gold, titanium)				✓	✓	✓	✓
M1,4 prosthetic screw on abutment (gold, titanium)	✓	✓					
M1,6 prosthetic screw (titanium)				✓			✓

All the values are indicative. In choosing the torque, follow the recommendations of the implant manufacturers.

VI. USING THE TORQUE WRENCH

The torque wrench is supplied unlubricated and unsterilised. Before first using it, the torque wrench must be lubricated and sterilised. (See "Hygiene and maintenance" §.).



Before use, check that the torque wrench is not damaged and no part of it is missing.

6.1. Connecting / disconnecting the rotary instrument




The use of protective gloves is recommended each time the tools are handled. Check

the condition of the rotating instruments to be used, and handle them properly and carefully.

Inspect the push button regularly to ensure that it is properly fixed.

To insert and tight the mandrel

- Apply pressure with the thumb to the center of the push button.
- Put the instrument into position by turning it slightly in both directions, until it engages, then release the push button.
- Check that the rotating instrument is properly attached, through a slight axial movement on each tool change.

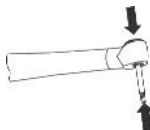


Removal of the rotating instrument

- Apply pressure with the thumb to the center of the push button.
- At the same time, remove the rotating instrument.

6.2. Torque selection


- Turn the selector A in the direction of the arrow while maintaining the B ring fixed.
- The torque values are indicated on the instrument.
- The implant manufacturer's recommendations are to be followed when choosing the torque.

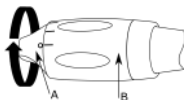

6.3. How to screw and reach the torque

Screwing/ Unscrewing

- Maintain the torque wrench in one hand.
- Turn the cylindrical component B until an audible click is heard, indicating that de-clutching has taken place and that the selected torque has been attained.

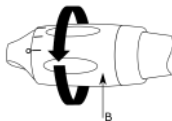

Medical device lifecycle

If used in a proper manner, all medical device parts have a two-year lifecycle (corresponding to 250 sterilisation cycles).

However, these indications are not a warranty because wear may appear prematurely, depending on how the device is maintained (cleaning and sterilisation).

VII. HYGIENE AND MAINTENANCE

The re-sterilization of reusable medical devices must be performed by a properly trained and protected member of staff, in accordance with the regulations currently in force.

The re-sterilization protocol must take into account the risk of infection.

Refer to the manufacturer's instructions for each product used. Take

particular care to adhere to the instructions regarding concentrations, lengths of exposure, the replacement of solutions, and the lifespan of the products. Never combine products. Adhere to the guidelines regarding the disposal of used products.

 Wear suitable protective clothing.

In order to avoid any risk of infection and injury, it is vital to wear protective gloves.

 Never immerse the torque wrench in any kind of solution.
Never clean the torque wrench in an ultrasonic tank.

7.1. Cleaning

The device must be cleaned immediately after the operation.

Use only detergent or disinfectant solutions with a neutral or slightly alkaline PH. It is extremely unadvisable to use products which may fix X proteins to the unit (eg.alcohol, aldehydes...).

 Do not use sodium hypochlorite (bleach) : there is a high risk of corrosion.

External cleaning

The torque wrench can be cleaned by using a brush under tap water. External disinfection using a spray or a wipe soaked with disinfectant products microbiologically controlled (ph between 2.5-9). In any case, it is important to follow the reaction time of the disinfectant product used.

 Do not use disinfectant containing alkaline and chlorine.
Do not immerse the torque wrench in a disinfectant solution or in an ultrasound container.

7.2. Lubrification

This should be done before using for the first time, and before each sterilisation.

 Keep away from any source of heat or ignition.
In particular, do not smoke!

Wear a protective face mask.

- Place a turbine-type nozzle (P/N 1933X) on the spray can.
- Place the spray can nozzle into the opening located on the backside of the torque wrench.
- Maintain the spray can nozzle against the torque wrench until the pressure drops again. Vaporise 1 to 2 seconds. The oil passes over the entire instrument, acting as a cleanser and a lubricant. Prior to each sterilisation, lubrication is to be carried out using a spray can.



7.3. Sterilisation



Read the instruction leaflet provided by the steriliser manufacturer. Adhere to the space specified between pouches and do not overload the steriliser.

- Only sterilise instruments which have previously been disinfected, cleaned, lubricated and tested.
- Remove the rotating instrument from the torque wrench prior to sterilisation. (See §6.1)
- Make sure that the device does not have any areas of corrosion or cracks, and check that it is operating properly. Ensure that the product is dry; if necessary, dry any residual water with medical quality pressurised air.
- Use sterilisation pouches suitable for the torque wrench and the steriliser.

Always use one single torque wrench per pouch.

- In order to prevent any retention of water, place the pouch in the steriliser in such a way that any concave parts are face down.

- If the steriliser has several types of cycle, choose a cycle designed for medical devices (a minimum of 135°C at 2.13 bars [275°F at 30.88 psi] for 18 minutes.)

- After each sterilization cycle, check that there is no water remaining on the inside and outside of the packaging. Make sure that the flow indicator has changed to the correct color.

- Keep the devices in the sterilisation pouches away from light, moisture and any kind of contamination. Follow the manufacturer's recommendations as seen on the packaging.

- The duration for which the device is kept after sterilisation should not exceed 1 month. Label the devices, specifying the expiration date.

After the expiration date, repeat the cleaning and sterilisation procedure.

VIII. CALIBRATION

In order to maintain the precision of the torque values, it is recommended that the torque wrench be returned for calibration every 2,000 uses or every 2 years.

IX. REPAIR

In the event of breakage, please contact your approved distributor or the Anthogyr after sales service department directly.



All repairs must be carried out with parts and subassemblies approved by the manufacturer.

Repairs must only be carried out by an approved distributor or by the after sales service department in the factory.

For service or repair, the equipment must be returned complete and sterile to the after sales service department. The equipment must be sent together with a document stating the problem as well as the complete name and address of the dentist.

For any claims to be considered under the guarantee, please enclose a copy of the invoice or the delivery note with the equipment.

For every service or repair, send the complete torque wrench.

Replacement parts are available for 7 years beyond sales discontinuation of the model.

X. GUARANTEES

This MD is guaranteed parts and labour against all manufacturing defects for 12 months from the date of invoice.

This guarantee does not apply to wear and tear parts.

All changes or additions to the product without the express agreement of Anthogyr render this guarantee null and void.

The guarantee becomes null and void if the technical instructions are not followed.

Anthogyr cannot be held responsible for damage resulting from or which could result from normal wear, use, cleaning or incorrect maintenance, the non-observance of instructions for use or connection, scaling or corrosion, impurities in the water supply system or unusual chemical or electrical influences or non observance of the instructions, maintenance instructions and assembly of Anthogyr and other manufacturer's instructions.

Delivery charges incurred when sending an instrument back to Anthogyr for repair will be paid by the client, even if the repair itself is covered by the guarantee.

Postage and packing fees when returning the instrument to the client are covered by the guarantee.

So that guarantee requests are taken into consideration, please attach a copy of the invoice or a copy of the delivery slip to the MD.

XI. MEDICAL DEVICE REFERENCES

Description

DYNAMOMETRIC C.A TORQCONTROL
DYNAMOMETRIC C.A TORQCONTROL
SPRAY ENDPIECE FOR TURBINE

References

15501
15501-P
1933X

XII. DISPOSAL OF THE PRODUCT

As far as is currently known, the product does not contain any substances which are harmful to the environment. The product must be sterilised before being disposed of. Observe national rules with regard to disposal.

Français

English (UK)

English (US)

Deutsch

Español

Italiano

Português

Nederlands

中國的

عربي



Anthogyr

2237 avenue André Lasquin
74700 SALLANCHES - FRANCE

Tél. +33 (0) 4 50 58 02 37

Fax +33 (0) 4 50 93 78 60

N°SAV / Repairs : 33 (0) 4 50 58 50 53

E-mail : contact@anthogyr.com

www.anthogyr.com

CE 0459